

Hamilton T1 Ventilator Setup for BiPAP/CPAP

Items Needed



T1 Ventilator




T1 Circuit



Non-vented Mask




Attach circuit and perform initial Pre-op tests

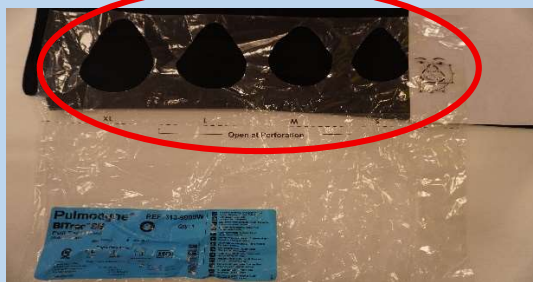


Determine appropriate sized mask for patient. Size guide on each mask package



Press Modes

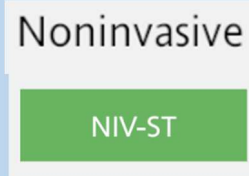

Select NIV-ST

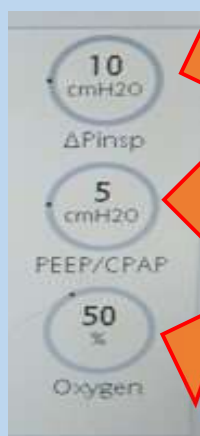

ΔPinsp

PEEP

FiO₂

REMEMBER the T1 ventilator is additive in pressure modes. The ΔPinsp is the **change** in pressure from PEEP to where you want the IPAP to be. Example: IPAP 15 PEEP 5, the ΔPinsp would be 10.

Setup Video Link



Setup

- Attach the circuit to the ventilator and perform the preop checks.
- Ensure that the elbow of the mask is NOT VENTED
- Determine the patient's current pressure settings (obtain from RT or copy from current floor device)
 - IPAP (if patient is BiPAP)
 - EPEP/PEEP/CPAP
 - FiO₂
- Press the **MODES** button
- Select **NIV-ST** mode of ventilation under the Noninvasive modes.
- Press **CONFIRM** to lock in the mode.
- Input the patient's current pressure settings previously obtained.
- Press **CONFIRM**.
- The ventilator will go back to the STANDBY screen and is ready for assisted ventilation.
- When ready, disconnect the patient from the hospital BiPAP/CPAP device and attach circuit to the mask.
- Secure the mask to the patient's face/head.
- Press START VENTILATION and monitor patient for changes through transport.
- Monitor and adjust for air leaks in mask. **Ensure proper seal.**