

# Hamilton T1 Ventilator Setup for Heated Humidified BiPAP/CPAP

## Items Needed



T1 Ventilator



Non-vented Mask



Sterile Water for Inhalation



Dual Limb Humidifier Circuit



Hamilton H900 Heater Humidifier



Hamilton Flow Sensor



Hamilton Expiratory Valve

## Heated Humidified BiPAP Therapy

Heated and humidified BiPAP/CPAP therapy is the most appropriate therapy when patients require prolonged noninvasive ventilation support. The main benefit of use is to improve compliance and comfort during the transport. It has also been shown to be effective and therapeutic for patients with thick secretions most associated with pneumonia or upper respiratory infection.

The principles of BiPAP/CPAP therapy on the Hamilton T1 Ventilator do not change, however there is additional equipment required to provide the heat and humidification. A dual limb circuit with a water chamber capable of attachment to the Hamilton H900 Heater is required. This circuit DOES NOT come prepackaged with a flow sensor or expiratory valve. Both are required to provide the therapy. Sterile water for inhalation is also required.

# Hamilton T1 Ventilator Setup for Heated Humidified BiPAP/CPAP




Limb connections.  
Blue on top White  
on bottom

Flow Sensor  
Connections.  
Blue to Blue,  
Clear to White



H900 Mode Selector

Connect Flow Sensor  
to the Y connection  
of the Circuit

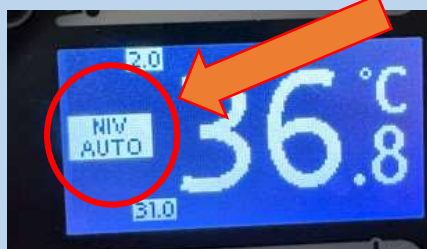

Setup Video  
Link



## Setup

- Attach the limb and flow sensor connections to the vent and dual limb circuit as show.
- Perform the preop checks.
  - Leak Test
  - Flow Sensor Test
- **Be sure not to remove the cap on the water tubing spike before performing the test. If removed, the pressure leak test will fail**
- Once the preop checks are performed, spike to sterile water bag and allow the chamber to fill to the “Max” line. (It will stop automatically when at that level.)
- Power on the H900 Heater device.
- Insert the water chamber into the H900 ensuring good contact with the sensors. If done correctly, two green arrows will appear.
- On the H900 device, choose the “NIV” mode. This will ensure proper heating and humidification levels.
- Allow the water chamber to heat to the preset temperature.
- Secure the H900 device to the EMS stretcher IV Pole.
- **The H900 MUST always be lower than the patient’s head!!**
- Input the BiPAP/CPAP settings in the same fashion as normal.
- Apply the mask to the patient initiate the therapy.
- Continuously monitor patient throughout transport.