

Hamilton T1 Ventilator Setup for Intubated Patients

Items Needed



T1 Ventilator



T1 Circuit



In-Line ETCO2 Sensor

Setup


- Attach the circuit to the ventilator and perform the preop checks.
- Obtain and enter the patient's height in CM (the T1 will automatically calculate the IBW).
- Determine the patient's current ventilation settings (obtain from RT or copy from the current floor ventilator).
 - Tidal Volume
 - PEEP
 - Respiratory Rate
 - FiO₂
- Press the **MODES** button.
- Select **(S)CMV+** mode of ventilation under the volume controlled modes of ventilation.
- Press **CONFIRM** to lock in the mode.
- Input the patient's current ventilation settings previously obtained.
- Press **CONFIRM**.
- The ventilator will go back to the STANDBY screen and is ready for ventilation.
- When ready, disconnect the patient from the hospital vent and attach circuit to the ETT.
- Press START VENTILATION and monitor patient for changes through transport.



Attach circuit and perform initial tests




Determine and input patient height

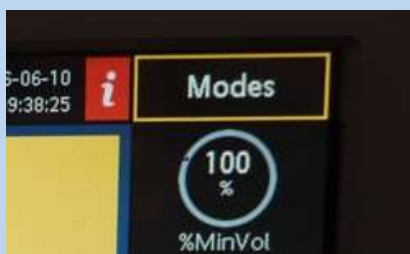

Press Modes



Select (S)CMV+

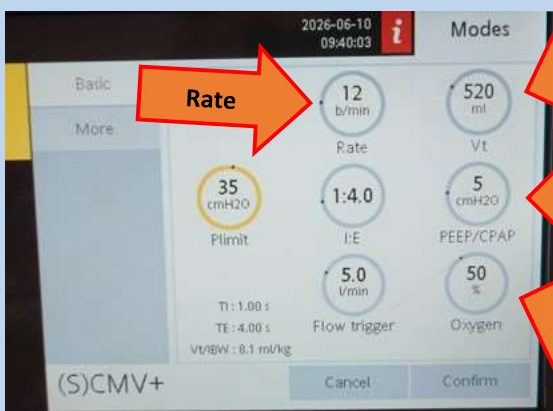

Video Setup Link

