


Nuvance Health EMS

Adult Heated Humidified High Flow NC



1. Choose appropriate sized cannula
 - should NOT occlude > 50% of nares
2. Connect Hamilton H900 inflow circuit to outflow port on vent.
3. Connect Optiflow cannula to the long outflow circuit tubing
4. On the Hamilton T1 choose “Non-invasive ‘high flow’” mode.
5. On the Hamilton H900 heater choose “High Flow” mode.
6. Set flow rate and FiO2.
7. Press “Start therapy”

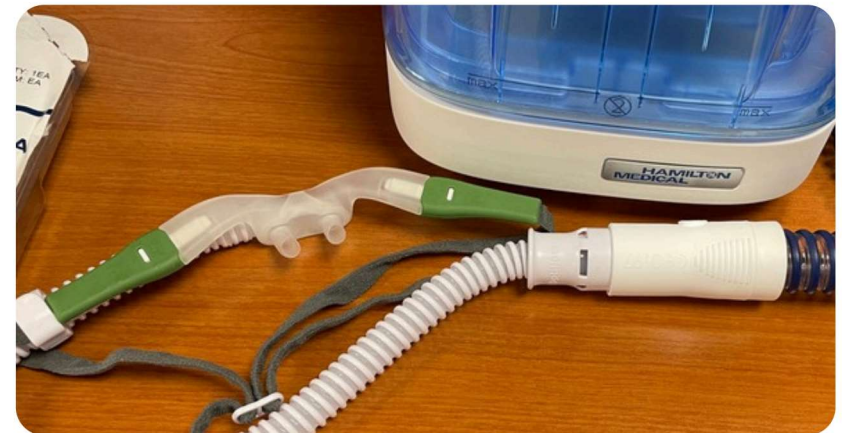

***Humidifier must be secured lower than the patient on stretcher ***

Flow rates

S 10-50 L/min


M 10-60 L/min

L 10-60 L/min




Setup Video
Link

2

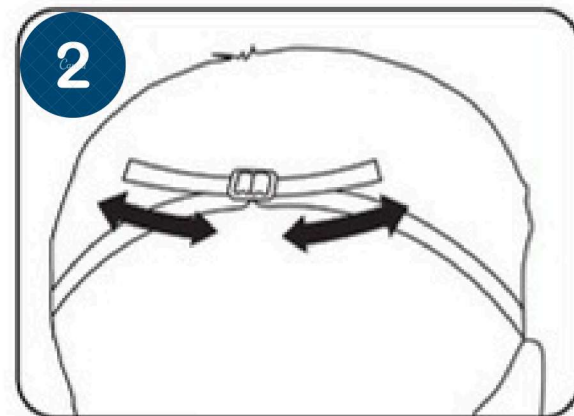

4

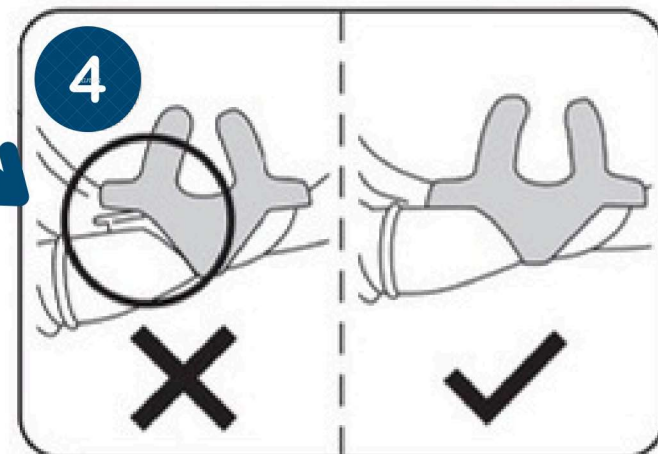

Connect to H900
outflow tubing


Benefits of HF

- Narrower tubes (compared with CPAP)
- Increase resistance to flow

- Delivers prescribed flow

- Interface designed to be unsealed

Optiflow+

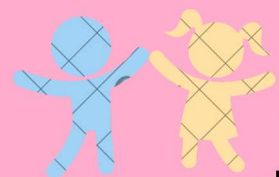


Nuvance Health EMS

Pediatric Heated Humidified High Flow NC



1. Choose appropriate sized cannula
 - should NOT occlude > 50% of nares
2. Connect Hamilton H900 inflow circuit to outflow port on vent.
3. Connect cannula to the long outflow circuit tubing using **Hudson RCI 15mm x 22mm adapter** and supplied adapter in the **Optiflow Jr 2 Ventilator Transition Kit**
4. On the Hamilton T1 choose “Non-invasive high flow” mode.
5. On the Hamilton H900 Heater choose “High Flow: Mode.
6. Set flow rate and FiO₂.
7. Press “Start therapy”

Setup
Video
Link

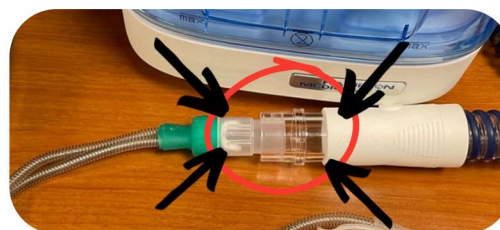


Humidifier must be secured lower than patient on stretcher



Fitting

- Select cannula size**
Select the appropriate cannula size; a clear gap must be visible around each prong.
- Prepare the skin**
Prepare the pt's skin according to protocol.
- Connect to gas source**
Connect the cannula to the gas source.
TIP: Place hand close to nasal prongs to ensure that there is gas
- Remove the first F&PWigglepad 2 tabs**
Remove the first backing tabs from the F&P Wigglepads 2 and avoid touching the adhesive.
- Insert cannula**
Insert the cannula into the nares. Ensure the cannula bridge rests close to the nose without touching the septum. **DO NOT** stretch the cannula during application. Stick the F&P Wigglepads 2 to the pt's cheeks.
- Secure cannula**
Remove the second backing tabs and stick the F&P Wigglepads 2 onto the cheeks.

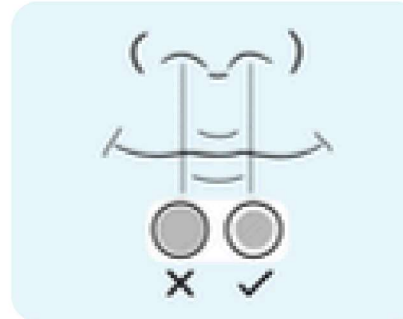

Ensure Gap Around Prongs



Flow rate

S		1.5-10 L/min
M		1.5-10 L/min
L		1.5-28 L/min
XL		1.5-31 L/min
XXL		1.5-36 L/min

Sizing

	S OJR412	
	Weight (kg)*	0.9 1 3.5 4
	Correlated age**	27 wkGA 28 wkGA 40 wkGA 42.5 wkGA
	M OJR414	
	Weight (kg)*	1 1.5 8 10
	Correlated age**	28 wkGA 31 wkGA 6.6 mo 15.4 mo
	L OJR416	
	Weight (kg)*	3 3.5 18 20
	Correlated age**	37.5 wkGA 40 wkGA 4.9 yr 5.6 yr
	XL OJR418	
	Weight (kg)*	5 7 25 30
	Correlated age**	47.5 wkGA 4.7 mo 7.6 yr 9.5 yr
	XXL OJR520	
	Weight (kg)*	8 12.5 30 40
	Correlated age**	6.4 mo 2.2 yr 9.5 yr 12 yr