

# Invasive Monitoring Setup Guide

## Items Needed



**Lifepak 35 with advanced monitoring capability**




**Lifepak 35 Invasive Pressure Monitoring Cables**



**Edwards Pressure Monitoring Set**



**Edwards Style Connector**

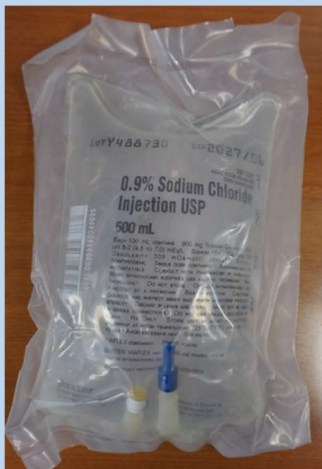

**500 mL Bag (Normal Saline)**



**500 mL Pressure Infuser**



**Velcro Armband Transducer Holer**

This guide will apply to any type of fluid filled invasive line monitoring. Types of monitoring include, but are not limited to:

- Arterial Pressure
- Central Venous Pressure
- Intracranial Pressure
- Pulmonary Artery Catheter (Swan-Ganz)
- Intra-Abdominal Pressure

AT THIS TIME NORTHWELL AMBULANCE OF CONNECTICUT IS ONLY AUTHORIZED TO TRANSPORT A PATIENT WHO REQUIRES ARTERIAL PRESSURE MONITORING. IF ANY OTHER TYPE OF MONITORING IS REQUIRED A HIGHER LEVEL OF CARE PROVIDER MUST ACCOMPANY THE EMS CREW

**Video Setup Link**



## Invasive Monitoring Setup Guide

### Setup

1. Confirm the sending facility's connection type. If current connection matches the Edwards style connector, skip to step # 18.
2. Spike the Saline bag and fill the drip chamber.
3. Place the bag into the pressure infuser and pump up the pressure to 300 mmhg.
4. Turn the 3-way stopcock attached to the transducer off toward the distal end of the line and pull the blue Snap-Tab to allow fluid to flush through the open cap.
5. Replace the open-ended clear cap with a solid white cap.
6. Turn the stopcock off to the solid cap.
7. Turn the distal in-line stopcock off to the end of line.
8. Pull the blue Snap-Tab to allow fluid to flush through the open cap.
9. Replace the open-ended clear cap with a solid white cap.
10. Turn the stopcock off to the solid cap.
11. Pull the blue Snap-Tab to allow fluid to flush through the distal end of the line.
12. Replace with the last remaining solid white cap.
13. Confirm that the pressure infuser remains at 300 mmhg. If needed pump up.




Edwards Style  
Connector




Snap-Tab

Select one

Stopcock

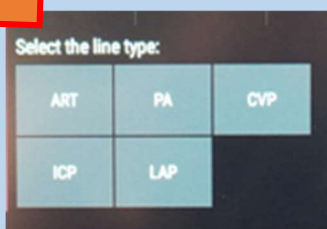

Press to Zero

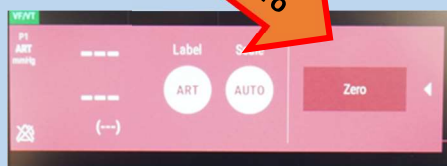

14. Locate the stopcock closest to the patient on the hospital's line set and close it toward the patient.
15. Remove the hospital's line set leaving the stopcock behind.
16. Attach the new line set to the stopcock ensuring that it is a wet-to-wet connection. **NO AIR PRESENT.**
17. Close the stopcock toward the cap.
18. Secure the new line transducer into the Dale armband holder.
19. Place the armband in line with the patient's phlebostatic axis and secure the Velcro.
20. Plug in the LifePak 35 invasive monitoring cables into the **RED** port on the monitor and connect the end of the transducer to a line.
21. Choose the appropriate monitoring desired.
22. Close the stopcock at the transducer toward the patient.
23. Remove the cap to allow the transducer to sample outside air pressure.
24. Press the channel on the monitor screen and go to the settings.
25. Press ZERO.
26. Once zeroing is complete, close the stopcock to the cap and confirm waveform and number readings.